

AMBER NA

Setting Up Your System



Quick Start Guide

To set up your system, follow the steps below. For more detailed information, see the AMBER NA User Guide.

Select a central location in your home for the AMBER NA where you spend most of your time. The AMBER NA can be placed on a table, desk, or counter. It also can be mounted on the wall.


After you complete steps 1 and 2, your AMBER NA will automatically call your Central station to guide you through the rest of the setup.

Note: Certain remote areas of your residence may lie outside the reception range of the AMBER NA, resulting in the unit not being able to receive transmissions directly. In this case, a Repeater (MCX-600 RL) can be used to increase the AMBER NA reception range.

Contact a Visonic representative for further details.


1. Connect the Base Unit to your Telephone Line:

A phone cord is already connected to the "wall jack" on your AMBER NA. Plug the other end of that phone cord into the phone jack on your wall.

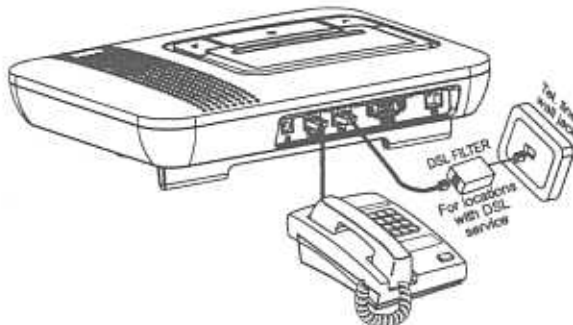


Note: Be aware of other phone line services such as DSL. If DSL service is present on the phone line, you must install a filter. For more information regarding the DSL filter, contact the service provider.

INSERT THE SAFETY CATCH TAB HORIZONTALLY AS SHOWN. TURN THE LEFT SIDE OF THE CATCH COUNTER-CLOCKWISE UNTIL IT FACES THE SCREW HOLE AND THEN FASTEN THE SAFETY CATCH WITH THE SCREW

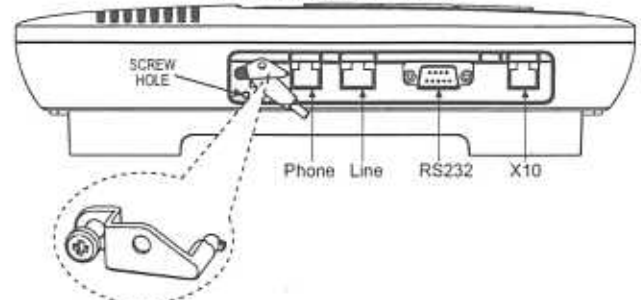

Note: The safety catch is optional, can be purchased upon request.

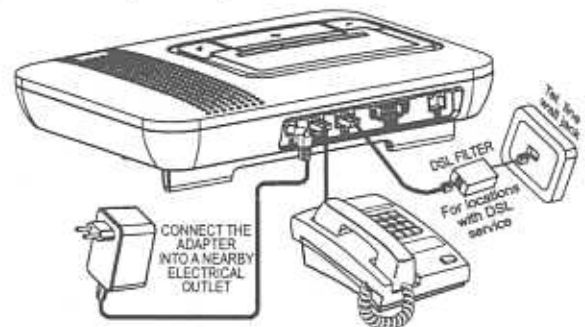

WARNING!

DO NOT USE AN OUTLET CONTROLLED BY A WALL SWITCH.

When electrical power is supplied to the Base Unit, the green light on the front panel illuminates.

2. The AMBER NA is now Connecting you to the Central Station:

After installation, a Monitoring Center representative will begin speaking to you.

When your Monitoring Center representative answers, you can hear the representative talking to you through the AMBER NA Unit. You can talk to the representative by speaking normally. You do not need to use the phone.

The Monitoring Center representative will guide you through the rest of the system setup.



WEEE Product Recycling Declaration

For information regarding the recycling of this product you must contact the company from which you originally purchased it. If you are discarding this product and not returning it for repair then you must ensure that it is returned as identified by your supplier. This product is not to be thrown away with everyday waste. Directive 2002/96/EC Waste Electrical and Electronic Equipment.



VISONIC LTD. (ISRAEL): P.O.B 22020 TEL-AVIV 61220 ISRAEL. PHONE: (972-3) 645-6789, FAX: (972-3) 645-6788

VISONIC INC. (U.S.A.): 65 WEST DUDLEY TOWN ROAD, BLOOMFIELD CT. 06002-1376. PHONE: (860) 243-0833, (800) 223-0020. FAX: (860) 242-8094

VISONIC LTD. (UK): FRASER ROAD, PRIORY BUSINESS PARK, BEDFORD MK44 3WH. PHONE: (0870) 730-0800 FAX: (0870) 730-0801

INTERNET: www.visonic.com

©VISONIC LTD. 2006 AMBER NA QUICK START GUIDE D-300495 (REV. 1, 09/06)

D-300495



MADE IN ISRAEL